


PAMS: Camera Quick Clamp

USER GUIDE

Introduction

The Camera Quick Clamp is an 3D printed adapter for attaching an assistive device to a section of 1/2 inch diameter PVC piping. The clamp has a 1/4" - 20 UNC camera thread stud that is compatible with various assistive devices and/or camera mount adapters. The clamp can easily be installed, adjusted, and removed by hand.

Features

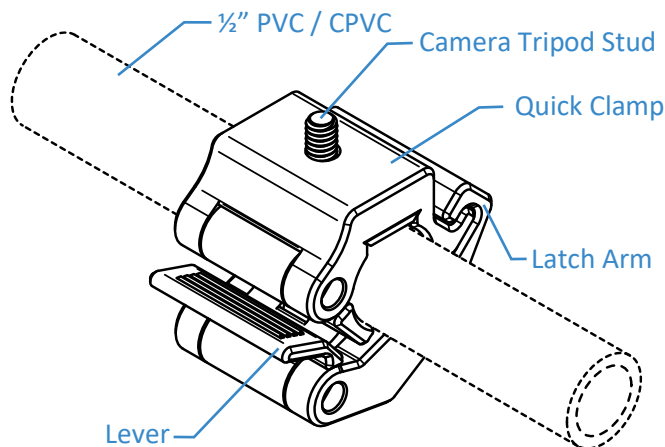


Specification


Envelope	Ø68 x 40 mm
Compatible Pipe	1/2 inch PVC 1/2 inch CPVC
Disinfection	Ethyl Alcohol, Isopropyl Alcohol, Hydrogen Peroxide

Usage


1. Form the quick clamp around the PVC piping.

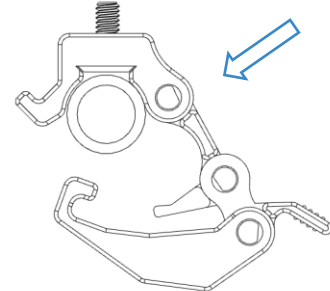

2. Position the latch arm onto the hook portion of the main body.

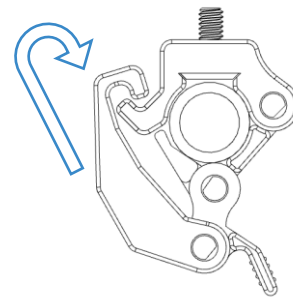

3. Press lever towards main body to tighten clamp on PVC.

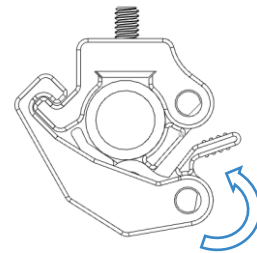

4. Attach assistive device using camera tripod stud.

